




Source of information: leaflet (note: the instructions might differ according to the brand of the acquired inhaler)

How to use the capsules with the inhaler

Do not swallow the capsules. You will breathe-in (inhale) the medicine in the capsules from the inhaler. Keep your inhaler dry. Handle with DRY hands. Do not use capsules with any other capsule inhaler, and do not use the inhaler to take any other capsule medicine.


Using the inhaler

1

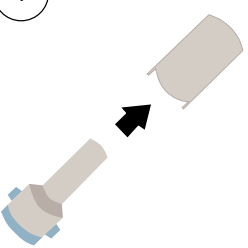


Do not remove a capsule from the blister card until you are ready for a dose. Pull off the inhaler cover.

2

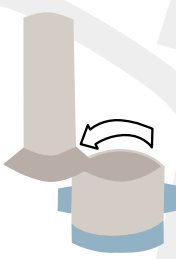


Hold the base of the inhaler firmly and twist the mouthpiece in the direction of the arrow to open. Push the buttons in on each side to make sure that you can see 4 pins in the capsule well of the inhaler.

3




Separate one capsule blister, peel the paper backing that covers one capsule and push it through the foil. Place the capsule in the capsule-chamber; never place it directly into the mouthpiece.

4




Twist the mouthpiece back to the closed position.

5

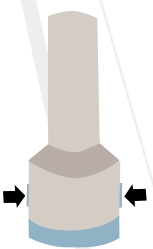


Hold the mouthpiece of the inhaler upright and press both buttons at the same time to pierce the capsule.

6

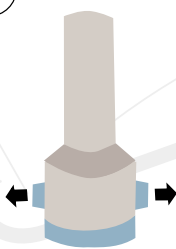


Release the buttons. Do not push the buttons a second time.

7



Breathe out (exhale) fully. Do not exhale into the inhaler mouthpiece.

8



Tilt your head back slightly. Keep the blue buttons to the left and right (not up and down). Place the mouthpiece in your mouth and close your lips around the mouthpiece.



9



Breathe in quickly and deeply. This will cause the capsule to spin around in the chamber and deliver your dose of medicine. You should hear a whirring noise.

10



Remove the inhaler from your mouth. Continue to hold your breath as long as you can.

11




Then exhale.

12


If you didn't hear any vibrating sound from the capsule oscillating in the chamber, you should review the steps in the leaflet that comes with the device and/or seek for help with a health professional.

13

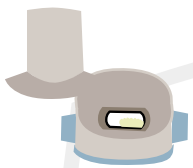


Open the inhaler to see if any powder is still in the capsule. If any powder remains in the capsule, you should review the steps in the leaflet that comes with the device and/or seek for help with a health professional.

14

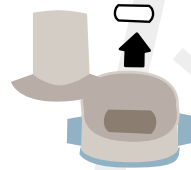

After use, open the inhaler, remove and discard the empty capsule. Do not leave a used capsule in the chamber.

15

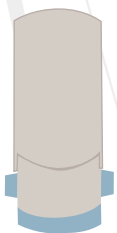

Close the mouthpiece and replace the cover.