




Source of information: leaflet (note: the instructions might differ according to the brand of the acquired inhaler)

### Testing your inhaler

If this is a new device or if you have not used it for 2 weeks or more, it must be tested before use by releasing 2 puffs into the air in the following way: 1) take the cover off from the mouthpiece, by pulling down on the lip at the back; 2) point the mouthpiece away from you so that the puffs of medicine will go into the air; 3) push the lever up so that it stays up. Then, to release a puff, push the dose release slide on the bottom of the device in the direction indicated by the arrow on it. To release the second puff you need to return the lever to its down position and follow steps 2 and 3 again. After releasing the second puff, return the lever to its down position, to be ready to take a puff of your medicine. Do not use the dose release slide to take your medicine; the Autohaler device will automatically release a dose when you begin to breathe in from the mouthpiece.


### Using your inhaler

1

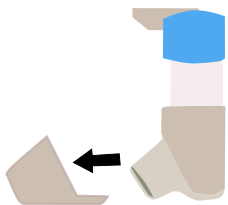


Take the cover off from the mouthpiece, by pulling down on the lip at the back.

2

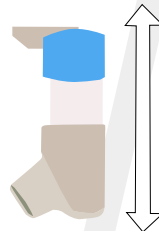

Shake the inhaler vigorously. Hold your Autohaler device upright as shown.

3

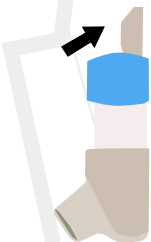

Push the lever up so that it stays up. Continue to hold your device upright, making sure that your hand is not blocking the air at the bottom of the device.

4



Breathe out as far as is comfortable.

5



Immediately place the mouthpiece in your mouth, and close your lips around it.

6



Breathe in slowly and deeply through the mouthpiece. Do not stop breathing in when you hear the slight click and feel the puff in your mouth as it is important that you carry on breathing in after the puff is released.




7



Hold your breath for 10 seconds.

8

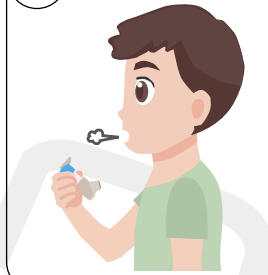

Breathe out slowly.

9



The lever must be lowered to the down position after each puff.

10



After use, replace the cover on the mouthpiece.

