




Source of information: leaflet (note: the instructions might differ according to the brand of the acquired inhaler)


Using your inhaler

1

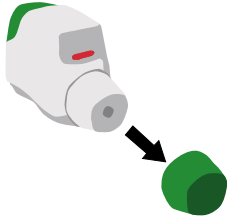


When you are about to take your dose of medicine, remove the protective cap by squeezing the arrows marked on each side and pulling outwards.

2

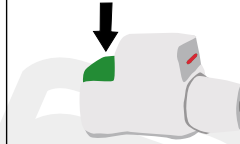


Before bringing the inhaler to your mouth, press the button all the way down.

3

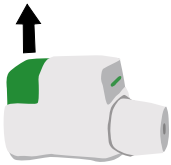

Then release it. **DO NOT CONTINUE TO HOLD THE BUTTON DOWN.**

4

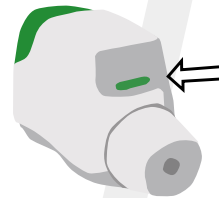

Make sure the coloured control window has changed to green. The green control window confirms that your medicine is ready for inhalation. **IF THE COLOURED CONTROL WINDOW STAYS RED, PLEASE REPEAT PRESS AND RELEASE ACTIONS.**

5



Before bringing the inhaler to your mouth, breathe out completely. Do not breathe out into the inhaler.

6



Put your lips tightly around the mouthpiece of the inhaler.

7



Breathe in **STRONGLY** and **DEEPLY** through your mouth. **ATTENTION: DO NOT HOLD THE GREEN BUTTON DOWN WHILE YOU ARE INHALING.** While you breathe in you will hear a "CLICK" which signals that you are using the inhaler correctly. Keep breathing in even after you have heard the "click" to be sure you get the full dose.



8




Remove the GENUAIR inhaler from your mouth and hold your breath for as long as is comfortable.

9



Breathe out gently.

10

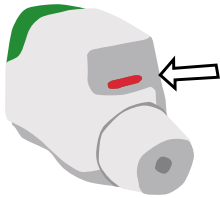

Make sure the control window has turned to red. This confirms that you have inhaled your full dose correctly. IF THE COLOURED CONTROL WINDOW IS STILL GREEN, PLEASE REPEAT INHALING STRONGLY AND DEEPLY THROUGH THE MOUTHPIECE.

11



Replace the protective cap by pressing it back onto the mouthpiece.