

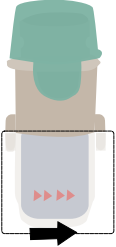


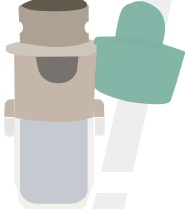
Source of information: leaflet (note: the instructions might differ according to the brand of the acquired inhaler)


Preparing your inhaler


Keep the cap closed. Press the safety catch while firmly pulling off the clear base with your other hand. Be careful not to touch the dose releasing button. Write the discard by date on the label (3 months from the date the cartridge is inserted). Insert the narrow end of the cartridge into the inhaler. Place the inhaler on a firm surface and push down firmly until it clicks into place. Put the clear base back into place until it clicks. Do not remove the clear base or the cartridge after it has been put together. Keep the cap closed. Turn the clear base in the direction of the arrows on the label until it clicks (half a turn). Open the cap until it snaps fully open. Point the inhaler towards the ground. Press the dose releasing button. If you do not see a mist, repeat steps from turning the clear base until a mist is seen. After a mist is seen, repeat steps from turning the clear base. After complete preparation of your inhaler, it will be ready to deliver the number of puffs on the label.


Using your inhaler


- 


Keep the cap closed. Turn the clear base in the direction of the arrows on the label until it clicks (half a turn).
- 


Open the cap until it snaps fully open.
- 

Breathe out slowly and fully.
- 

Close your lips around the mouthpiece without covering the air vents. Point the inhaler to the back of your throat.
- 

While taking a slow, deep breath through your mouth, Press release button and continue to breathe in.
- 

Hold your breath for 10 seconds or for as long as comfortable.
- 

Breathe out gently.
- 

Close the cap until you use your inhaler again.